



Bands



Global Experts in Orthodontics

03

GC Orthodontics
GLOBAL EXPERTS IN ORTHODONTICS

Bands Summary

Band coding	103
Sure-Snap molar bands	106
A-Fit molar bands	107
Sure-Snap bicuspid bands	108





Bands

We are offering two types of high quality molar bands. Both are manufactured in Japan by Tomy® Inc

Band coding

Sure-Snap molar bands How are P/Ns built?




UL18



SUL26




A-Fit molar bands How are P/Ns built?

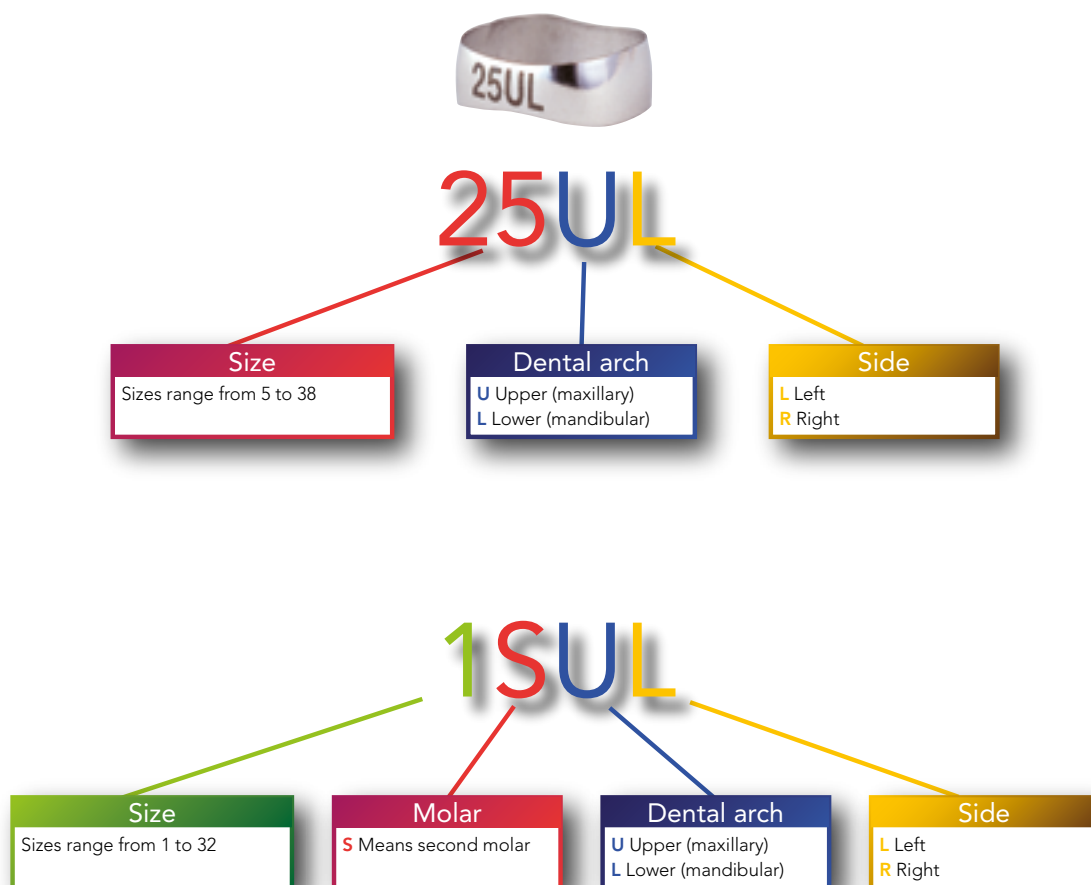



Sure-Snap bicuspid bands How are P/Ns built?

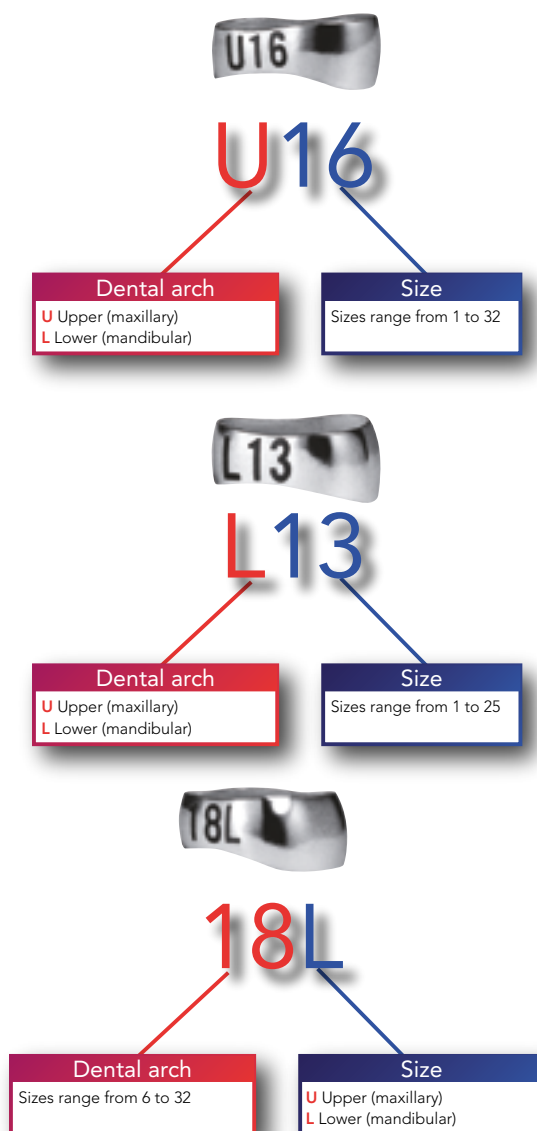



Sure-Snap molar bands

The alloy Sure-Snap is made for makes it both hard and elastic. The difference in thickness between the occlusal edge and the gingival edge makes it:

- Resistant to the pressure of the band-pusher (occlusal edge)
- Elastic enough to get over the contour of the teeth without any deformation (gingival edge).

With Sure-Snap, you will get a special anatomy for each type of molar, an easy-to-read permanent marking and an electroeroded intrados for unique retention.

They are available in 32 sizes for 1st molars; and in 30 sizes for 2nd molars

Sure-Snap Kits first molars - Distribution by quadrant: (UR, UL, LR, LL)

Size	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32
Starter											1	2	2	2	2	2	2	2	2	2	2	2	2									
Intro	1	1	1	1	1	1	1	1	1	1	2	2	2	2	2	3	3	3	3	2	2	2	2	2	1	1	1	1	1	1	1	
Standard	1	1	1	1	2	2	2	3	4	4	4	5	5	6	6	6	6	6	6	5	5	4	3	2	2	2	1	1	1	1	1	
Master	1	1	1	2	2	2	4	4	6	6	8	10	10	12	14	16	16	16	14	12	10	8	6	4	4	2	2	2	2	1	1	1

Sure-Snap Kits second molars - Distribution by quadrant: (UR, UL, LR, LL)

Size	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32
Starter											1	2	2	2	2	2	2	2	2	2	2	2	2									
Intro	1	1	1	1	1	1	1	1	1	2	2	2	2	2	2	3	3	3	3	2	2	2	2	2	2	1	1	1	1	1		
Standard	1	1	1	2	2	2	2	3	4	4	4	5	5	6	6	6	6	6	6	5	5	4	3	2	2	2	2	1	1	1		
Master	1	1	1	2	2	2	4	4	6	6	8	10	10	12	14	16	16	16	16	12	10	8	6	4	4	2	2	2	2	1		

Upper first molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 32 Sizes

Standard Kit: 200 Assorted bands, 32 Sizes

Master Kit: 400 Assorted bands, 32 Sizes

Upper second molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 30 Sizes

Standard Kit: 200 Assorted bands, 30 Sizes

Master Kit: 400 Assorted bands, 30 Sizes

Lower first molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 32 Sizes

Standard Kit: 200 Assorted bands, 32 Sizes

Master Kit: 400 Assorted bands, 32 Sizes

Lower second molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 30 Sizes

Standard Kit: 200 Assorted bands, 30

Sizes Master Kit: 400 Assorted bands, 30 Sizes



A-Fit molar bands

The A-Fit Band is featuring all the qualities of Sure-Snap. The difference is a higher crown height for the A-Fit Band and a larger choice of sizes in the higher numbers. A-Fit Bands identification: see page 109.

A-Fit Kits first molars - Distribution by quadrant: (UR, UL, LR, LL)

Size	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32	33	34	35	36	37	38	
Starter									1	2	2	2	2	2	2	2	2	2	2	2	2														
Intro	1	1	1	2	2	2	2	2	2	2	2	3	3	3	2	2	2	2	2	1	1	1	1	1	1	1	1	1	1	1	1	1			
Std.	1	2	2	3	4	4	4	5	6	6	6	6	6	6	5	5	4	4	3	3	2	2	2	2	1	1	1	1	1	1	1				
Master	2	4	4	6	8	8	8	10	12	12	12	12	12	12	10	10	8	8	6	6	4	4	4	4	2	2	2	2	2	2	2				

A-Fit Kits second molars - Distribution by quadrant: (UR, UL, LR, LL)

Size	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32	
Starter												1	2	2	2	2	2	2	2	2	2	2	2										
Intro	1	1	1	1	1	1	1	1	1	1	1	2	2	2	2	2	2	3	3	3	3	2	2	2	2	1	1	1	1	1	1	1	1
Standard	1	1	1	1	2	2	2	3	4	4	4	5	5	6	6	6	6	6	6	5	5	4	3	2	2	2	1	1	1	1	1	1	1
Master	1	1	1	2	2	2	4	4	6	6	8	10	10	12	14	16	16	16	14	12	10	8	6	4	4	2	2	2	2	2	1	1	1

Upper first molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 34 Sizes

Standard Kit: 200 Assorted bands, 34 Sizes

Master Kit: 400 Assorted bands, 34 Sizes

Upper second molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 32 Sizes

Standard Kit: 200 Assorted bands, 32 Sizes

Master Kit: 400 Assorted bands, 32 Sizes

Lower first molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 34 Sizes

Standard Kit: 200 Assorted bands, 34 Sizes

Master Kit: 400 Assorted bands, 34 Sizes

Lower second molars, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 32 Sizes

Standard Kit: 200 Assorted bands, 32 Sizes

Master Kit: 400 Assorted bands, 32 Sizes



Sure-Snap bicuspid bands

Manufacturing quality is identical for these bands and for our molar bands.

Band setting is made through a rocking movement, alternating mesial and distal pressures. They are available in three different anatomies:

- 1st maxillary bicuspid (32 Sizes)
- 1st mandibular bicuspid (25 Sizes)
- 2nd mandibular bicuspid (27 Sizes)

Thanks to their easy setting, you will not hesitate any more to use them instead of brackets for severe rotations or in case of bond failures.

Sure-Snap Kits upper bicuspid - Distribution

Size	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32
Starter						2	2	2	2	2	2	2	2	2	2	4	4	4	2	2	2	2	2	2	2	2						
Intro	2	2	2	2	2	2	2	2	2	2	4	4	4	4	6	6	6	6	6	4	4	4	4	2	2	2	2	2	2	2	2	2
Standard	2	2	2	2	2	4	6	6	6	6	6	8	8	10	12	14	14	14	12	10	8	8	6	6	6	4	4	4	2	2	2	2
Master	4	4	4	6	6	10	10	14	14	16	16	16	18	18	20	22	22	22	20	18	16	16	14	14	12	12	8	8	8	4	4	4

Sure-Snap Kits 1st lower bicuspid - Distribution

Size	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	
Starter			2	4	4	4	4	4	4	4	4	4	4	4	4											
Intro	2	2	2	2	2	2	4	4	4	4	4	6	6	8	8	8	8	4	4	4	4	2	2	2	2	
Standard	2	2	2	2	4	4	8	8	8	8	10	12	12	18	18	18	18	10	8	8	8	4	4	2	2	
Master	4	4	4	6	6	8	8	14	14	18	18	24	24	36	36	36	36	24	24	16	16	8	8	4	4	

Sure-Snap Kits 2nd lower bicuspid - Distribution

Size	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32
Starter						1	1	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	1	1
Intro						2	2	2	2	2	2	4	4	4	6	6	8	8	8	8	6	6	4	4	4	2	2	2	2	2	2	2
Standard						2	2	2	4	4	6	6	6	8	10	12	14	16	16	16	14	12	10	8	6	6	6	4	4	2	2	2
Master						2	4	4	6	8	10	12	12	16	20	26	28	34	34	34	28	26	20	16	14	12	10	8	6	4	4	2

Upper bicuspid, bands only

Starter Kit: 50 Assorted bands, 22 Sizes

Intro Kit: 100 Assorted bands, 32 Sizes

Standard Kit: 200 Assorted bands, 32 Sizes

Master Kit: 400 Assorted bands, 32 Sizes

1st lower bicuspid, bands only

Starter Kit: 50 Assorted bands, 13 Sizes

Intro Kit: 100 Assorted bands, 25 Sizes

Standard Kit: 200 Assorted bands, 25 Sizes

Master Kit: 400 Assorted bands, 25 Sizes

2nd lower bicuspid, bands only

Starter Kit: 50 Assorted bands, 27 Sizes

Intro Kit: 100 Assorted bands, 27 Sizes

Standard Kit: 200 Assorted bands, 27 Sizes

Master Kit: 400 Assorted bands, 27 Sizes